

Polycarbonate Vial Sets

For high throughput processing of seeds, leaves and animal tissues

WISBIOMED

120 Boothbay Ave.
San Mateo CA, 94404


TEL (650) 242-8368

FAX (650) 242-8244

Email: info@wisbiomed.com


Sample Preparation Products

24 Well PC Vial Set (Optional Lined Caps)

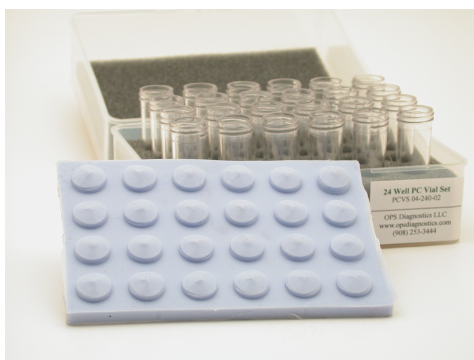

Designed to the SBS 24 well format, the 24 Well Vial Sets are for high throughput homogenization of kernels, seeds and animal tissues. Each set contains 24 polycarbonate screw cap tubes, 24 stainless steel 3/8 inch grinding balls, and 24 polypropylene caps. Vials can be sealed with caps or an optional sealing mat (contact us for details) during grinding and then capped and handled as individual samples. Tubes are formatted correctly for liquid handling. They are also available with lined caps for tightly sealing vials. Vials with lined caps have been used effectively used to contain viruses during tissue homogenization.

15 mL PC Vial Set

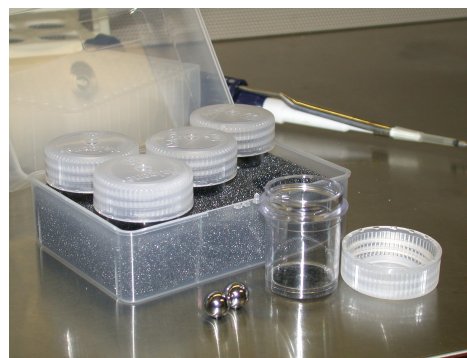

The 15 ml Vial Set was developed for researchers who required pooled samples for analysis. This product was designed for the difficult job of dry grinding 15 lyophilized corn kernels for field trials. Each 15 ml Vial Set comes with 5 polycarbonate vials per box with two stainless steel grinding balls per vial.

Product Number	Description	Quantity	Cost
OP-PCVS 15-050-01	15 ml PC Vial Sets	Case of 10	\$249
OP-PCVS 04-240-02	24 Well PC Vial Sets	Case of 10	\$249
OP-PCVS 04-240-03	24 Well PC Vial Sets with lined Caps	Case of 10	\$270

Related Products

OP-PEVS 04-240-04	24 Well PE Vial Sets	Case of 10	\$202
OP-PEVS 04-240-05	24 Well PE Vial Sets with lined Caps	Case of 10	\$223

WISBIOMED

120 Boothbay Ave.
San Mateo CA 94404