

EVERLAST™ ROCKER 247


Benchmark 



Model BR5000

- + Everlast™ 5 year warranty
- + Stainless steel surface
- + 16 lbs. load capacity
- + Adjustable speed & tilt angle
- + Optional stacking platform

**Extra large
platform, 16x12 in.
(40 x 33 cm)**

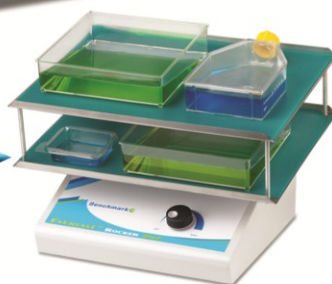

Model BR5000 shown with BR5000-STACK

Technical Data:

Speed Range:	3 to 80 rpm
Tilt Angle:	0 to 20° (user adjustable)
Platform Size:	16 x 12 in. / 40.5 x 30cm
Load Capacity:	16 lb / 7.5 kg
Dimensions:	16 x 13 x 6 in. 40.5 x 33 x 15.25 cm
Operating Temp.:	0°C to +65°C
Weight:	20 lb / 9 kg
Electrical:	115V or 230V, 50-60Hz
Warranty:	5 years

With an unprecedented 5 year limited warranty, Benchmark's Everlast 247 Rocker is built to last. Twenty four hours a day, seven days a week... No problem! The powerful motor and heavy-duty drive system remain reliable and consistent, even under constant usage in a cold room (min. 0°C) or incubator (max. 65°C).

The rocking speed and tilt angle can easily be adjusted to vary between gentle and aggressive agitation. Both the standard platform and the (optional) stacking platform are constructed of stainless steel, preventing against any rusting or corrosion. The Everlast 247 is designed to provide years of uninterrupted service.

Ordering Information:

BR5000* Everlast™ Rocker 247 with non-slip rubber mat

Optional Accessories:

BR5000-STACK Optional stacking platform for BR5000
Clearance: 3.5 in.

* To order a product in 230V / 50Hz, please add -E to item number.

Distributed by WISBIOMED LLC
1-866-692-1249 wisbiomed.com